

SHIN-NIPPON
by **Rexxam**

CHART PROJECTOR

CP-500

CHART PROJECTOR CP-500



50 charts and various masks

Landolt, Letter, Child, Snellen, 50 kinds of charts, and Horizontal line, Vertical line, one chart, various masks can be selected by the remote control box.

Fast, silent chart disk rotation

It moves to the following chart in a maximum of 0.3 seconds.

User Program mode

2 kinds of programs can be selected. Each program can store 48 kinds (maximum) of chart and mask.

Various Switches

- 1 Lamp switch
- 2 No use
- 3 Program switch (1,2)
- 4 Chart/Mask switch
- 5 Horizontal line switch
- 6 Vertical line switch
- 7 One chart switch
- 8 Back switch
- 9 Next switch

SPECIFICATIONS


Number of charts and masks	50 charts and 38 masks
Luminance	min.80-max.320Cd/m2(200Cd/m2 at 5m)
Chart and mask selection speed	0.3seconds
Program	2 kinds(48 charts or masks)
Tilting angle	±10°
Auto power off	After 5 min.
Lamp	White LED
Projection distance	2~7 m
Projection magnification	30 x (at 5m)
Projection size	335mm x 222mm (at 5m)
Dimensions	W248xL336xH234mm
Weight	4.5kg
Power source	100~240V
Power consumption	30VA

STANDARD ACCESSORIES

REMOTE CONTROL BOX, SUM4 BATTERY X 2, SCREEN, TEST SCALE SHEET, DUST COVER, SPARE FUSE

OPTIONAL ACCESSORIES

TABLE BASE STAND, WALL MOUNT, CHART PROJECTOR ARM, POLA FILTER GLASS, R&G GLASS, DEVIATOR

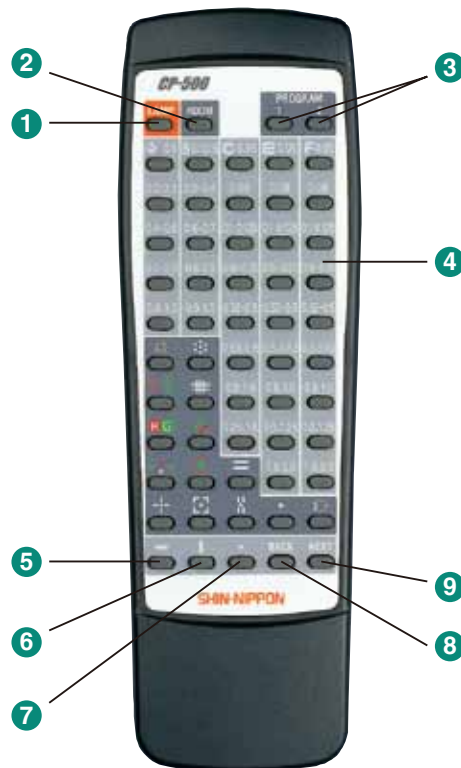


Adjustment of the brightness

The brightness of the projector can be changeable by the control switch according to the brightness of the operation room.



A type chart (B type chart is also available.)

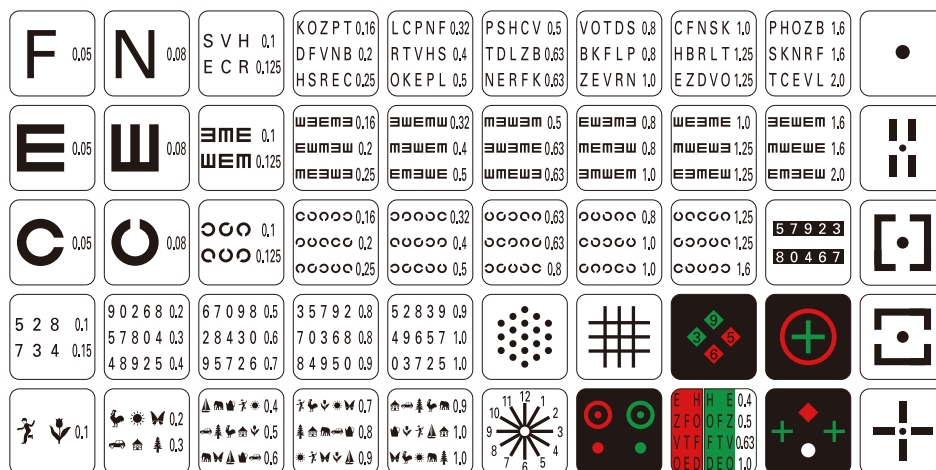

Design and specifications are subject to change without notice.